

This Page Is Inserted by IFW Operations  
and is not a part of the Official Record

## **BEST AVAILABLE IMAGES**

Defective images within this document are accurate representations of the original documents submitted by the applicant.

Defects in the images may include (but are not limited to):

- BLACK BORDERS
- TEXT CUT OFF AT TOP, BOTTOM OR SIDES
- FADED TEXT
- ILLEGIBLE TEXT
- SKEWED/SLANTED IMAGES
- COLORED PHOTOS
- BLACK OR VERY BLACK AND WHITE DARK PHOTOS
- GRAY SCALE DOCUMENTS


IMAGES ARE BEST AVAILABLE COPY.

**As rescanning documents *will not* correct images,  
please do not report the images to the  
Image Problems Mailbox.**



**PCT**WORLD INTELLECTUAL PROPERTY ORGANIZATION  
International Bureau

## INTERNATIONAL APPLICATION PUBLISHED UNDER THE PATENT COOPERATION TREATY (PCT)
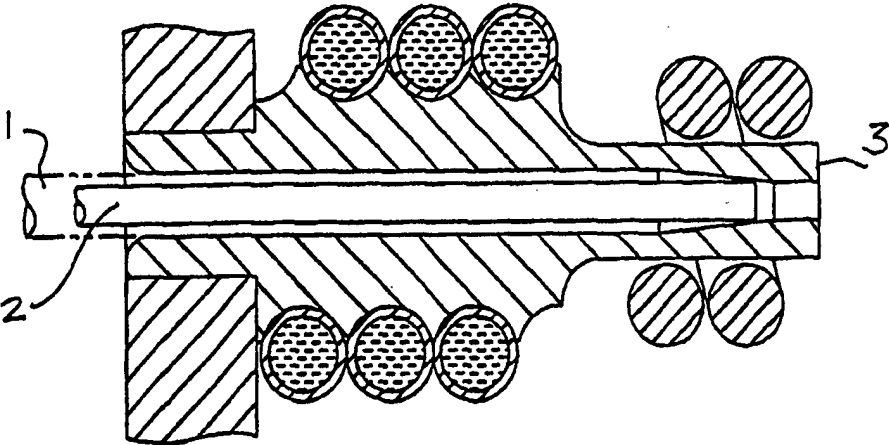
<b>(51) International Patent Classification <sup>6</sup> :</b> <b>A61L 29/00</b>	<b>A1</b>	<b>(11) International Publication Number:</b> <b>WO 99/47187</b> <b>(43) International Publication Date:</b> 23 September 1999 (23.09.99)
<b>(21) International Application Number:</b> PCT/US99/05846 <b>(22) International Filing Date:</b> 17 March 1999 (17.03.99)  <b>(30) Priority Data:</b> 09/040,758 18 March 1998 (18.03.98) US  <b>(71) Applicant:</b> BECTON, DICKINSON AND COMPANY [US/US]; 1 Becton Drive, Franklin Lakes, NJ 07417-1880 (US).  <b>(72) Inventors:</b> KHAN, Azhar, J.; 3031 Maple Way, West Valley City, UT 84119 (US). KHAN, Mohammad, A.; 10005 South Falcon Hurst Drive, Sandy, UT 84092 (US).  <b>(74) Agent:</b> LEE, Eric, M.; Becton, Dickinson and Company, 1 Becton Drive, Franklin Lakes, NJ 07417-1880 (US).		<b>(81) Designated States:</b> AU, BR, CA, CN, IN, JP, MX, SG, European patent (AT, BE, CH, CY, DE, DK, ES, FI, FR, GB, GR, IE, IT, LU, MC, NL, PT, SE).  <b>Published</b> <i>With international search report. Before the expiration of the time limit for amending the claims and to be republished in the event of the receipt of amendments.</i>
<b>(54) Title:</b> WATER SOLUBLE LUBRICANT FOR MEDICAL DEVICES    <b>(57) Abstract</b> <p>This invention relates to a new water soluble lubricant for a medical device, such as a catheter and an introducer needle. The lubricant is a silicone based water soluble compound. Preferably, a silicone surfactant is used such as a block copolymer polyalkylene oxide-modified polydimethylsiloxane. Alternatively, an alkylene oxide modified silicone glycol may be used as the lubricant. In addition, lecithin can be added as a lubricant. The lubrication solution includes a solution stabilizer to clarify the solution and antimicrobial agents to inhibit microbial growth in the water solution or on the coated product. Vitamin E or its derivatives may also be used in the lubrication solution.</p>		

**FOR THE PURPOSES OF INFORMATION ONLY**

Codes used to identify States party to the PCT on the front pages of pamphlets publishing international applications under the PCT.

AL	Albania	ES	Spain	LS	Lesotho	SI	Slovenia
AM	Armenia	FI	Finland	LT	Lithuania	SK	Slovakia
AT	Austria	FR	France	LU	Luxembourg	SN	Senegal
AU	Australia	GA	Gabon	LV	Latvia	SZ	Swaziland
AZ	Azerbaijan	GB	United Kingdom	MC	Monaco	TD	Chad
BA	Bosnia and Herzegovina	GE	Georgia	MD	Republic of Moldova	TG	Togo
BB	Barbados	GH	Ghana	MG	Madagascar	TJ	Tajikistan
BE	Belgium	GN	Guinea	MK	The former Yugoslav Republic of Macedonia	TM	Turkmenistan
BF	Burkina Faso	GR	Greece	ML	Mali	TR	Turkey
BG	Bulgaria	HU	Hungary	MN	Mongolia	TT	Trinidad and Tobago
BJ	Benin	IE	Ireland	MR	Mauritania	UA	Ukraine
BR	Brazil	IL	Israel	MW	Malawi	UG	Uganda
BY	Belarus	IS	Iceland	MX	Mexico	US	United States of America
CA	Canada	IT	Italy	NE	Niger	UZ	Uzbekistan
CF	Central African Republic	JP	Japan	NL	Netherlands	VN	Viet Nam
CG	Congo	KE	Kenya	NO	Norway	YU	Yugoslavia
CH	Switzerland	KG	Kyrgyzstan	NZ	New Zealand	ZW	Zimbabwe
CI	Côte d'Ivoire	KP	Democratic People's Republic of Korea	PL	Poland		
CM	Cameroon	KR	Republic of Korea	PT	Portugal		
CN	China	KZ	Kazakstan	RO	Romania		
CU	Cuba	LC	Saint Lucia	RU	Russian Federation		
CZ	Czech Republic	LI	Liechtenstein	SD	Sudan		
DE	Germany	LK	Sri Lanka	SE	Sweden		
DK	Denmark	LR	Liberia	SG	Singapore		
EE	Estonia						

5

10

### **WATER SOLUBLE LUBRICANT FOR MEDICAL DEVICES**

This application is a continuation-in-part of application Serial No. 08/798,810 filed February 12, 1997, which is a continuation of application Serial No. 08/670,296 filed June 27, 1996, which is a continuation of application Serial No. 08/294,212 filed August 22, 1994.

#### **Background of the Invention**

This invention relates in general to a novel lubricant used in connection with medical devices such as intravenous (IV) catheters.

IV catheters are designed to infuse normal IV solutions, including antibiotics and other drugs, into a patient. The typical catheter is hollow and is extruded out of a suitable plastic material such as Teflon, polyvinyl chloride, polyethylene, polyurethane or polyether urethane. These catheters are also used to withdraw blood from the patient for normal blood-gas analysis as well as other blood work. In order to insert an IV catheter into a patient, an introducer needle is used. The needle is typically stainless steel and is hollow. Its distal tip is ground to a sharp tip for easy insertion into the patient. The catheter is also hollow and is initially located coaxially around the introducer needle in an "over the needle" arrangement. The internal diameter of the catheter tip is slightly less than the outer diameter of the tip of the needle so that the catheter tip has an interference fit on the needle tip. The interference fit is necessary so that when the catheter and introducer needle assembly is taken out of the package, the catheter remains snugly on the needle and does not easily slip off. This interference fit also facilitates insertion of the catheter and introducer needle assembly

5 into the patient's vein because it minimizes the chance that the catheter tip will fold over or peel back on the needle tip.

Placement of the catheter and introducer needle into the patient may cause sharp pain to the patient. In order to facilitate insertion of the catheter and introducer needle into the vein and to minimize patient discomfort, the catheter and introducer  
10 needle can both be lubricated. Most IV catheters are lubricated with polydimethylsiloxane silicone fluid. However, some IV catheters are not lubricated at all.

The polydimethylsiloxane silicone fluid may be applied to the surface of the catheter and needle by wiping the surfaces with the lubricant. Alternatively, the  
15 catheter and needle can be separately dipped into a solution of polydimethylsiloxane silicone fluid and a solvent. This is generally the preferred method of applying the lubricant because a consistent, controlled and uniform coating can be achieved. The polydimethylsiloxane silicone fluid must be dissolved in an organic solvent because the silicone oil in this compound is hydrophobic. Typically, the solution contains 2.5%  
20 silicone oil. The catheter and needle are then separately dipped into this solution.

Until recently the preferred solvent has been freon because it is non-flammable and flashes off, i.e. evaporates, readily after the polydimethylsiloxane silicone fluid solution has been applied to the catheter and needle. Although freon has been preferred, it does suffer some drawbacks. For instance, the high evaporation rate of  
25 freon causes the polydimethylsiloxane silicone fluid to concentrate on the surface of the solution in which the catheter and needle are dipped. Because of this high evaporation rate, the solution is difficult to control. Moreover this solution is expensive because of the large loss of freon during the coating process. In addition, freon is a chlorofluorocarbon (CFC) which is thought to react with and destroy the earth's  
30 protective ozone layer. Thus the manufacture and use of such CFC's will eventually cease in the near future.

5           Other solvents will have to be used in order for silicone fluid to be applied to catheters and introducer needles as a lubricant. Other solvents include alcohol and hydrocarbons. However, alcohol and hydrocarbons are highly combustible and are therefore too dangerous for use in manufacturing.

10           Previous efforts have been made to provide a lubrication system that uses water as the solvent. Such systems typically are a three part system that uses a lubricant, a surfactant and water. U.S. Patent No. 5,266,359 discloses such a three part system that uses water as the solvent. Although, such three part lubrication systems that use water as the solvent are "environmentally friendly" and are not flammable, they could be improved. Since a three part lubrication system by definition requires three separate  
15 elements, a certain amount of cost is associated with the materials and manufacturing processing needed to apply the system to a medical device to be lubricated.

          Another important issue concerning IV catheters is the shape of the catheter tip. The shape of the catheter tip must produce minimal trauma to the patient during insertion of the catheter into the patient and while the catheter is in place in the patient.  
20 Such a preferred tip shape that provides these characteristics has a tapered outer wall and an angled tip and is disclosed in U.S. Patent No. 4,588,398. A process for making that catheter tip is disclosed in U.S. Patent No. 4,661,300. In this process, the catheter is placed on a mandrel. A die having an interior molding surface, which is tapered according to the tip desired on the catheter, is aligned axially with the mandrel. The  
25 catheter tip is heated, typically using RF energy, so that it is flowable. The mandrel and die are brought together so the distal edge of the mandrel engages the tapered portion of the die. This action cleanly forms a smooth and uniform tapered tip for the catheter.

          After the catheter is tipped, it must be free of defects such as incomplete formation, substantial flash or jagged edges. The tip must also look smooth and be free  
30 of roll-overs. In addition, the length of the catheter must remain within a desired specification after the tipping process. If, during the tipping process, the thermoplastic

5 material sticks to the die or the mandrel, the length will vary greatly due to stretching and the tip will not be free of defects. Visual or microscopic examination may be used to determine if there are any tip defects and if the length of the catheter is within specifications.

Typically a lubricant is used to allow the tipped catheter to be easily removed  
10 from the mandrel and die. If a lubricant is not used, the tipped catheter could stick to the mandrel or die resulting in a deformed catheter when it is removed from the mandrel or die. Standard tipping lubricants include polydimethyl siloxanes such as Dow Corning DC 360 or curable silicones such as Dow Corning 44159 MDX which are amine terminated and moisture curable. Non-curable amine terminated polydimethyl  
15 siloxanes have also been used for this purpose. Such lubricants are described in, for example, U.S. Patent Nos. 3,574,673; 4,904,433; and 5,185,006. Again, these silicone oils used are problematic as discussed above.

#### **Summary of the Invention**

20 It is therefore an object of this invention to provide a lubrication solution for a medical device, such as a catheter and an introducer needle assembly, that is inexpensive and easy to control.

It is another object of this invention to provide a lubricant for a medical device, such as a catheter and an introducer needle assembly, that does not require the use of  
25 a CFC as a solvent.

It is yet another object of this invention to provide a lubrication solution for a medical device, such as a catheter and an introducer needle assembly, that is "environmentally friendly."

It is still another object of this invention to provide a lubrication solution for a  
30 medical device, such as a catheter and an introducer needle assembly, that is not flammable.

5 It is a further object of this invention to provide a lubrication solution that requires less material and cost than other lubrication solutions.

The lubrication solution of this invention is a two part system that uses water as the solvent. The lubricant is a silicone based water soluble material. No other material is needed for lubricity. Instead the silicone in the material alone provides the lubricity  
10 and the material is soluble in water. Thus, the solution of the silicone based water soluble material and water allows the silicone to be easily applied to the medical device or other apparatus and, when the water evaporates, ensures that the surface on which the material is applied is lubricious. Preferably, the silicone based water soluble material is a silicone surfactant or a silicone glycol. The silicone material is preferably  
15 non-ionic because it may have lower toxicity than the ionic form. The lubrication solution may also include vitamin E or its derivatives. In addition, the lubrication solution in which the device to be lubricated is dipped may include a solution stabilizer to clarify the solution and antimicrobial agents to inhibit microbial growth in the water solution or on the coated product.

20 The tipping process of this invention requires that the tip portion of the catheter blank be dipped in the solution or to otherwise have this solution applied to the surfaces to be lubricated such as by brushing or spraying the solution onto such surfaces. Once the water evaporates, the catheter blank is mounted on a mandrel and heated. A die and the mandrel are brought into engagement to form the catheter tip. The tipped  
25 catheter is then easily and quickly removed from the die and the mandrel.

The above and other objects and advantages of the invention will be apparent upon consideration of the following detailed description.

#### **Brief Description of the Drawings**

30 The preferred embodiments are illustrated in the drawings in which like reference numerals refer to like elements and in which:

5           FIG 1 is a cross-sectional view of a portion of a die and mandrel arrangement with a catheter blank thereon with the mandrel engaged with the die to form the catheter tip.


#### Detailed Description of the Invention

10           Although this invention is discussed in terms of its application to IV catheters and introducer needles, it is to be understood that this invention could be used on other medical devices where a lubricious surface on the device is desirable.

          This invention uses a silicone based water soluble compound as the lubricant since the silicone in the compound is a good lubricating fluid. The silicone based water  
15 soluble compound should comprise between about 0.25% to about 40.0%, preferably from about 2.0% to about 6.0% of the lubrication solution. The lubrication solution may also include vitamin E or its derivatives. Preferably between about 0.1% to about 1.0% of vitamin E or its derivative is used. The lubrication solution should also include a solution stabilizer to clarify the solution and antimicrobial agents to inhibit microbial  
20 growth in the water solution or on the coated product. The solution stabilizer should comprise between about 0.1% to about 10%, preferably between about 0.2% to about 1.0% of the solution. The antimicrobial agent should comprise between about 0.001% to about 5.0%, preferably between about 0.002% to about 0.05% of the solution.


          The lubrication solution of this invention is a two part system of the silicone  
25 based water soluble compound and water. One silicone based water soluble compound that can be used is a silicone surfactant such as a Silwet silicone surfactant. Silwet is the trade name of a class of silicone surfactants sold by OSI Specialities, Inc. These surfactants are polyalkylene oxide-modified polydimethylsiloxane block copolymers. They are similar to standard silicone fluids except the  
30 polydimethylsiloxane backbone has polyalkylene oxide side chains similar to non-ionic surfactants like poly-(oxyethylene) poly-(oxypropylene) block copolymers known as pluronic polyols. The side chains are terminated with either hydroxy or low alkoxy end

- 5 modified silicone glycols sold by PPG. These compounds have the following structural formula:

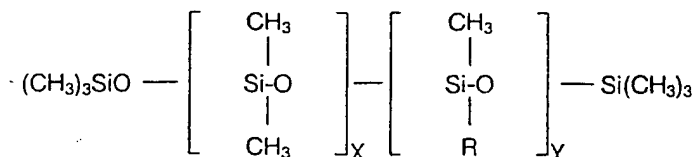

10


where R is defined as:

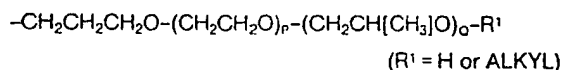

- MASIL® 1066D Silicone Surfactant is similar to the model above, but the silicone-organic linkage is through an Si-O-C bond, susceptible to hydrolysis under acidic or mildly basic conditions. The alkylene oxide modified silicone glycol can be used alone as the lubricant or it can be used in combination with the silicone surfactant.


- Vitamin E and its derivative, Vitamin E acetate, may also be added to the solution. Vitamin E, which is chemically known as alpha-tocopherol, is an antioxidant and thus prevents degradation of the solution. Vitamin E and its derivative, vitamin E acetate, are both oily products and also enhance the lubricity of this lubrication system.
- The molecular structure of vitamin E is given below:

25

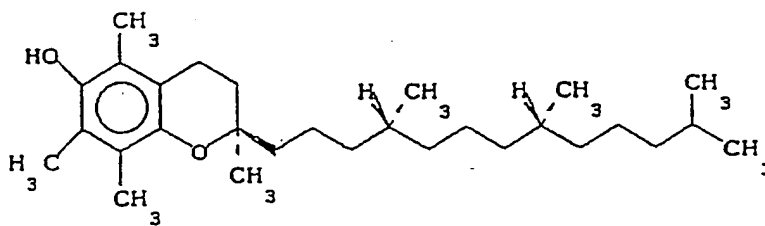

- Since the lubricant solution is prepared in water, it is highly desirable that a small amount of an effective antimicrobial agent be present to serve as a preservative. In the absence of such an agent, micro-organisms may grow in the solution and make the solution toxic. There are several commercial antimicrobial agents available. These are iodophors; phenols; phenolic compounds such as para-chloro-meta-xyleneol; biguanides

5 groups. One of these surfactants, Silwet L7001 has a molecular weight of 20,000 and a viscosity of 1700 centistokes. Its chemical formula is shown below:

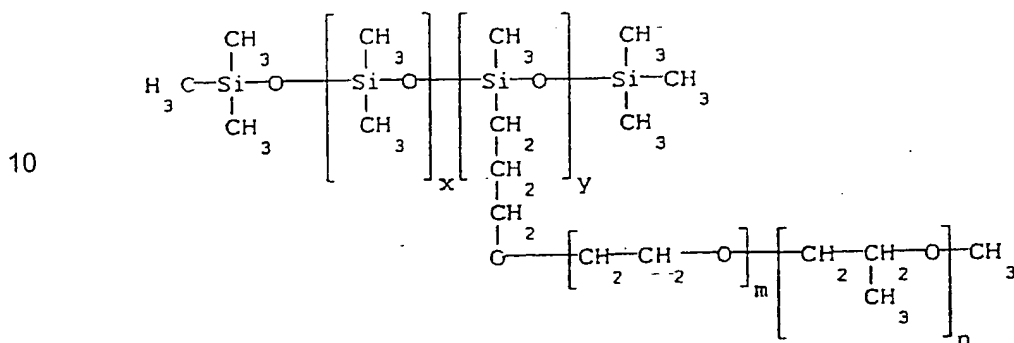

15 The amino-modified silicone polyether copolymer can also be used as the lubricant alone or in combination with the polyalkylene oxide-modified polydimethylsiloxane block copolymer. These surfactants are soluble in a wide variety of solvents such as CFC, alcohol, acetone, and water. The polyalkylene oxide chains also promote wetting on polyether urethane surfaces. Polyether urethane is a material that is used  
20 for making IV catheters. The similarity in chemical structure between these side chains and the soft segment of polyether urethane promote the surfactant's affinity for the catheter surface.

Although surfactants can be irritating or toxic depending on exposure levels, the Silwet silicone surfactants are copolymers of two polymeric materials, silicone fluids and polyalkylene oxides, which are low in toxicity. In particular, Silwet L7001 has a very low order of acute toxicity by swallowing, or skin penetration and is minimally irritating to the skin and is not irritating to the eyes. Thus, there should be no toxicity problems when the lubricant of this invention is used on a patient.

Another silicone based water soluble compound that can be used as the  
30 lubricant of this invention is an alkylene oxide modified silicone glycol. A preferred  
compound is sold under the trade name Masil which is a class of alkylene oxide

- 5           The following examples show the beneficial result of using the lubricant of this invention:

Example No. 1

- 10           Initial studies were conducted by separately dipping 20 gauge (ga) catheters and introducer needles into 4%, 8% and 32% Silwet L7001 silicone surfactant and water solutions and then assembling the catheters and needles. The assemblies were penetration tested through 13.5 mil. thick natural latex film.

15		<u>4% Silwet L7001</u>	<u>8% Silwet L7001</u>	<u>32% Silwet L7001</u>
	Needle tip (g)	24.5 (0.9)	23.3 (2.9)	24.3 (5.8)
	Needle transition (g)	18.5 (2.7)	17.0 (2.2)	15.0 (1.0)
	Needle heel (g)	11.0 (0.7)	10.5 (1.5)	8.0 (0.0)
	Catheter tip (g)	26.0 (3.7)	28.3 (5.3)	25.3 (1.9)
20	Catheter taper (g)	15.3 (0.8)	15.0 (3.0)	10.3 (0.8)
	Catheter drag (g)	7.3 (1.9)	4.5 (0.5)	4.0 (0.0)

NOTE: ( ) = standard deviation. Sample size = 4.


- 25           The values given represent the resistance in grams when the device is penetrated through the latex membrane. These values are comparable to those of currently marketed products that are lubricated with dimethylsiloxane fluids and are better than products that are unlubricated. This example thus shows that a silicone surfactant works as a lubricant for a catheter and/or introducer needle.
- 30

Example No. 2

- 20 gauge catheters and needles were separately dipped in a solution containing 8% Silwet L7001 silicone surfactant, 0.25% vitamin E, 0.5% vitamin E acetate, and 0.026% cosmocil. The catheters and needles were then assembled and these assemblies were penetration tested through 13.5 mil. thick natural latex film. The results are shown below:
- 35

5 such as chlorhexidine gluconate and polyhexamethylene biguanide hydrochloride (cosmocil). Cosmocil is used because it is less toxic than the other antimicrobial agents and is used as a preservative in contact lens cleaning solutions. The molecular structure of cosmocil is given below:

10

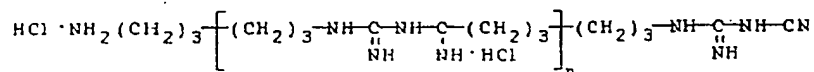

15

The molecular weight of this product is  $2100 \pm 300$ .

The process of forming the catheter tip will now be described. First, the lubricant solution is applied to the surface to be lubricated, i.e. the catheter blank. The catheter blank 1 to be tipped is then mounted onto a mandrel 2 supported on a carriage (not shown). A die 3 having an interior molding surface, with at least one portion of which is tapered according to the tip desired on the finished catheter, is aligned axially with the mandrel 2. The carriage is moved toward the die 3 such that the end of the catheter blank 1 to be tipped engages the interior molding surface. The carriage is halted at a point after the catheter blank 1 has engaged the tapered portions of the interior molding surface of the die 3, and is biased toward the die 3 with only sufficient force to cause the catheter blank 1 to move further into the die as the catheter blank 1 is heated by RF energy to its melting point and begins to flow. See FIG. 1. The catheter blank 1 is then allowed to heat in the die 3 and is moved into the die 3 as it begins to melt and flow under the biasing force. The mandrel 2 is positioned so that its distal corner is against the tapered portion of the die 3, thus cleanly forming the leading edge of the catheter. The die 3 and catheter therein are then cooled, and the carriage is reversed such that the catheter may be withdrawn from the die 3. Finally, the catheter is removed from the mandrel 2. There is no flash left on the leading edge because of the contact between mandrel 2 and die 3.

5

	<u>No Quaternary Ammonium Salt</u>	<u>1% Benzalkonium Chloride</u>
Needle tip (g)	24.3 (4.6)	24.0 (3.6)
Needle transition (g)	19.0 (1.0)	19.3 (1.4)
Needle heel (g)	11.3 (1.2)	10.0 (0.5)
10 Catheter tip (g)	19.8 (2.1)	24.2 (2.3)
Catheter taper (g)	15.5 (0.9)	17.3 (1.8)
Catheter drag (g)	4.8 (0.5)	4.9 (0.4)

NOTE: ( ) = standard deviation.

- 15 These results indicate that there is no adverse effect on the lubricant caused by using a quaternary ammonium salt in the lubrication system.

#### Example No. 4

- 20 Since amino-modified silicones are lubricious when applied to metal surfaces, the following lubricant solutions were prepared.

	<u>Lubricant I</u>	<u>Lubricant II</u>
Silicone Surfactant	4.50	4.50
25 Amino-modified silicone polyether copolymer	0.00	0.50
Vitamin E (%)	0.25	0.25
Cosmocil (PPM)	50	50
Water Q.S. to	100	100

30

20 gauge catheter products were assembled after dipping the catheter and the introducer needle separately into the above solutions. The tip adhesion force between the catheter and the needle was measured after aging the assemblies at 90°C for 2 1/2 days. The results are tabulated below:

35

	<u>Tip Adhesion (lbs)</u>
Lubricant I	0.348 (0.11)
Lubricant II	0.188 (0.03)

- 40 NOTE: ( ) = standard deviation.

5	Needle tip (g)	18.3 (1.9)
	Needle transition (g)	14.7 (0.6)
	Needle heel (g)	7.0 (0.3)
	Catheter tip (g)	20.2 (2.3)
	Catheter taper (g)	10.5 (0.6)
10	Catheter drag (g)	3.5 (1.3)

NOTE: ( ) = standard deviation. Sample size = 5.

It can be seen that the use of vitamin E and/or its derivatives as part of the lubricant in the solution improves the penetration forces.

15        When vitamin E and/or its derivatives are included in the silicone surfactant solution, a rather cloudy solution is obtained indicating that the solution is not homogenous. If a quaternary ammonium salt is added at a certain concentration, a clear solution is obtained. For example, Sylguard, sold by Dow Corning, which is a reactive quat and benzalkonium chloride can be added to the solution. Other  
20        quaternary ammonium salts, such as benzethonium chloride, could also be used.

### Example No. 3

      A solution containing 6.0% Silwet L7001 silicone surfactant, 0.36% vitamin E and  
25        1.0% benzalkonium chloride was used as the lubricant solution. In this example, the 1.0% benzalkonium chloride was in the form of 2.0% Hyamine 3500. Hyamine 3500 is a trade name for the benzalkonium chloride solution sold by Rohm and Hass and contains 50% of the active ingredient. As a comparison, a lubricant solution with no benzalkonium chloride was also used. 18 gauge catheters and needles were  
30        separately dipped and then assembled. These catheter assemblies were penetrated through dental dam (natural latex film).

	<u>Composition A</u>	<u>Composition B</u>	<u>Composition C</u>
Needle tip (g)	237.0 (63.7)	226.0 (68.6)	209.0 (46.1)
Catheter tip (g)	299.0 (53.5)	379.0 (42.7)	378.0 (41.2)
Catheter taper (g)	320.0 (37.9)	284.0 (34.0)	327.9 (39.4)
Catheter drag (g)	51.0(8.0)	46.0 (9.7)	42.0 (11.2)

NOTE: ( ) = standard deviation. Sample size = 5.

These results on thick latex membrane show that the amount of silicone surfactant applied to the catheter and the needle has opposite effects on catheter drag and needle tip penetration. The higher the silicone surfactant concentration, the lower the drag, but at high lubricant concentrations the needle tip penetration value increases. On fresh cow skin, the needle tip penetration values decreased with increased concentration.

#### Example No. 6

The following formulations were tested.

<u>Catheter Lubricant</u>		<u>Needle Lubricant</u>	
Silicone Surfactant	4.75% $\pm$ 0.25%	Silicone Surfactant	2.38% $\pm$ 0.25%
Amino-modified silicone polyether copolymer	0.525% $\pm$ 0.025%	Amino-modified silicone polyether copolymer	0.525% $\pm$ 0.025%
Vitamin E	0.263% $\pm$ 0.013%	Vitamin E	0.263% $\pm$ 0.013%
Cosmocil	50 ppm	Cosmocil	50 ppm
Water Q.S. to	100	Water Q.S. to	100

20 gauge catheters and needles were assembled as before and penetration tested through natural latex film of 13.5 mils thickness. At the same time, commercial 20 gauge Insyte® catheter products were tested under the same conditions for comparison. The results are tabulated below:

	<u>Experimental Sample</u>	<u>Commercial Sample</u>
Needle tip (g)	22.6 (8.6)	19.5 (5.7)
Needle transition (g)	17.3 (1.9)	13.4 (1.2)
Needle heel (g)	10.9 (1.1)	7.0 (0.6)
Catheter tip (g)	18.1 (3.0)	16.3 (2.4)
Catheter taper (g)	16.2 (2.1)	10.6 (1.0)
Catheter drag (g)	4.8 (1.9)	3.0 (0.6)

NOTE: ( ) = standard deviation. Sample size 10.

- 5 These test results indicate that the use of an amino-modified silicone-polyether copolymer does lubricate the metal surface better than the unmodified silicone surfactant alone.

10 Example No. 5

The following statistical experiment was run to optimize the formulation of the lubrication system of this invention. The following table lists some of the composition variations used.

	<u>Composition A</u>	<u>Composition B</u>	<u>Composition C</u>
15 Silicone Surfactant	3.0	4.5	6.0
Amino-modified silicone polyether copolymer	0.5	0.25	0.5
Vitamin E (%)	0.125	0.13	0.125
20 Cosmocil (PPM)	50	50	50
Water Q.S. to	100	100	100

- 20 gauge catheter assemblies were assembled after dipping the catheters and the introducer needle separately. After aging the assemblies at 90°C for 2 1/2 days, these
- 25 assemblies were penetration tested through thick natural latex film of 86.7 mils thickness. The thick latex was used because it is closer to the thickness of mammalian skin and because small differences in penetration forces are magnified. The results are tabulated below:

	<u>Composition A</u>	<u>Composition B</u>	<u>Composition C</u>
30 Needle tip (g)	156.0 (11.7)	150.4 (13.1)	155.9 (9.5)
Catheter tip (g)	90.6 (10.3)	87.0 (10.5)	98.3 (7.0)
Catheter drag (g)	31.6 (2.9)	32.0 (2.3)	27.6 (1.8)

- 35 NOTE: ( ) = standard deviation. Sample size = 8.

The above assemblies were also aged at 90°C for two weeks to simulate a five year shelf-life and penetration tested through fresh cow skin. The results are tabulated below:

- 5 systems of the invention was less than the clotting resulting from the commercial lubricant. The calculated risk of clotting was three times less for the water based lubricant versus the commercial lubricant.

### Example No. 9

- 10 Studies were conducted by dipping 18 gauge catheters and introducer needle assemblies into the various lubrication solutions shown below. The amounts of each ingredient is in percent by weight unless otherwise specified.

Ingredients	1	2	3	4	5	6	7	8
Alkylene oxide modified silicone glycol (Masil 1066D)	4.5	--	--	--	5.0	--	--	--
Alkylene oxide modified silicone glycol (Masil 1066C)	--	4.5	--	--	--	5.0	--	--
Alkylene oxide modified silicone glycol (Masil 29)	--	--	4.5	--	--	--	5.0	--
Alkylene oxide modified silicone glycol (Masil 280)	--	--	--	4.5	--	--	--	5.0
Amino modified silicone surfactant (Y12686)	0.5	0.5	0.5	0.5	--	--	--	--
Vitamin E	0.25	0.25	0.25	0.25	0.25	0.25	0.25	0.25
Cosmocil	50 ppm	50 ppm	50 ppm	50 ppm	50 ppm	50 ppm	50 ppm	50 ppm
Water	94.75	94.75	94.75	94.75	94.75	94.75	94.75	94.75

15

The assemblies were penetration tested through a latex film. The results are tabulated below:

	Comp. 1	Comp. 2	Comp. 3	Comp. 4	Comp. 5 <sup>1</sup>	Comp. 6	Comp. 7	Comp. 8
Needle tip (g)	21.95(3.2)	22.80(4.3)	20.45(1.8)	20.65(3.1)	--	22.95(4.2)	20.78(4.0)	22.70(2.5)
Catheter tip (g)	13.95(1.5)	12.65(0.7)	12.70(0.7)	13.85(1.0)	--	13.10(0.9)	14.15(1.1)	15.40(1.1)
Catheter drag (g)	3.30(0.4)	3.86(1.5)	3.90(1.2)	5.15(1.0)	--	4.15(1.0)	5.95(1.0)	7.44(1.2)

20

NOTE: ( ) = standard deviation. Sample size = 10.  
Compound 5 was not tested because of difficulties in mixing the solution.

- 5 As evident, the results are comparable for the commercial product and the experimental sample. The commercial sample used polydimethylsiloxane silicone fluid as the lubricant.

10 Example No. 7

- 22 gauge catheters and needles were assembled as before using the formulations described in Example 6 and penetration tested through sheep skin on the hind legs just below the knee. The penetration forces were measured using sensitive force transducers. The data were collected on a computer for analysis. For comparison, commercial 22 gauge catheter and needle assemblies using polydimethylsiloxane silicone fluid as the lubricant were used. The results are summarized below:

20		<u>Experimental Sample</u>	<u>Commercial Sample</u>
	Max. Catheter/Needle Tip (g)	142.57	164.67
	Max. Catheter drag (g)	70.31	271.68

- 25 As is evident, the new lubrication system is superior to the commercial product.

Example No. 8

- Polyurethane tubes were lubricated using the formulations of this invention described in Example No. 6 as well as polydimethylsiloxane silicone fluid. These tubes were implanted into the aorta of rabbits to determine the level of associated clotting and emboli formulations. After three days the animals were sacrificed and the tubes were examined while they were still in the aorta. Clotting, if present, was photographed and recorded as to size, number and location. The kidneys were also examined for renal infarcts which would indicate that emboli had formed and traveled downstream to lodge in the small arteries of the kidney. The clotting resulting from the water based lubricant

- 5 sticky in nature and will bond to the surface of a metal under compression. Since the catheter tip sits over the stainless steel needle with an interference fit, it will stick to the needle unless the needle is lubricated.

#### Example No. 11

- 10 The following lubricants were used during the tipping process as a tipping lubricant and to lubricate the needle and the catheter. The control group is included for comparison.

		<u>Tipping Lube</u>	<u>Needle Lube</u>	<u>Catheter Lube</u>
15	Silicone surfactant (%)	$2.375 \pm 0.25$	$2.375 \pm 0.25$	$4.75 \pm 0.25$
	Amino-modified silicone surfactant (%)	$0.525 \pm 0.025$	$0.525 \pm 0.025$	$0.525 \pm 0.025$
	Vitamin E (%)	$0.2625 \pm 0.0125$	$0.2625 \pm 0.0125$	$0.2625 \pm 0.0125$
20	Cosmocil (PPM)	50	50	50
	Water (%)	$96.8375 \pm 0.2875$	$98.8375 \pm 0.2875$	$94.4625 \pm 0.2875$
<b>Control</b>				
25	DC 12600 cstk Silicone (%)	-----	-----	2.0
	Masil 1 MM cstk Silicone (%)	-----	2.4	-----
	PS-513 Amino-modified Silicone (%)	0.5	-----	-----
30	HCFC-141b, Genesolve 2000 (%)	99.5	97.6	98.0

- 20 gauge catheter products were assembled by first using the catheter tipping lubricant for tipping the catheter, lubricating the needle and catheter separately using the respective lubricants and finally assembling the catheter assembly. The products were aged at 90°C for two weeks and tested for tip adhesion. The products were also penetration tested through latex membrane 13.5 mils thick. The results are described below:

- 5 These results are comparable to the results tabulated in Example No. 1.

The tip adhesion force between the catheter and the needle was measured for certain solutions and tabulated below.

Tip Adhesion	Comp. 3	Comp. 4	Comp. 6
Tip Adhesion (g) (non aged)	185.93 (23.37)	254.64 (25.20)	234.10 (23.40)
Tip Adhesion (g) (aged 60°C 18h)	132.19 (45.89)	209.99 (104.26)	428.70 (183.29)

10

NOTE: ( ) = standard deviation. Sample size = 10.

#### Example No. 10

The following table identifies the effects of using different amounts of a silicone

- 15 surfactant and an amino-modified silicone surfactant as the lubricant in a catheter tipping operation.

	<u>1</u>	<u>2</u>	<u>3</u>	<u>4</u>
20 Silicone surfactant (%)	3.00	3.00	4.50	3.00
Amino-modified silicone surfactant(%)	0.00	0.25	0.00	0.50
Vitamin E (%)	0.125	0.00	0.25	0.125
Cosmocil (PPM)	50	50	50	50
Water (%)	96.88	96.75	95.25	96.38
25 No. of Tips Sticking (out of 60)	8	0	1	0
Catheter Length (inches) (Target length is 2.0 inches)	1.99866	2.00030	1.99953	1.99993

- 30 If the catheter tip sticks, the catheter length will not be 2.0 inches due to stretching of the catheter tip as it is removed from the mandrel. As can be seen, when the lubricant provided adequate lubricity, no tips stuck to the mandrel and the catheter length was very close to the target length of 2 inches.

- 35 If the adhesion between the catheter tip and the introducer needle is too high, the needle cannot be easily withdrawn from the catheter after the assembly is placed into the patient's vein. Many thermoplastic materials such as polyurethanes are very

- 5 was held stationary at certain angle rotations to make sure that the needle is rubbing the wedge. The needle was pulled out and the resistive force was measured. The results are given below:

	<u>Product Tested</u>	<u>Resistive Force (lbs)</u>
10	Unlubricated 20(ga) catheter	0.400
	20(ga) Catheter Lubricated with Silicone	0.160
	20(ga) Catheter Lubricated with Lubricant of this invention as described above	0.078

- The data clearly shows the lubricant effectiveness of the water soluble lubrication system of this invention.

#### Example No.13

- A quaternary ammonium salt was used as a solution stabilizer for the tipping lubricant solution to ensure that the solution is homogenous. A lubricant containing 3% Silwet L7001, 0.5% Sylguard, which is a reactive quaternary compound, 0.25% Urea, and 96.25% water was used for tipping catheters. In all products tested, the tip quality was acceptable. Other quaternary ammonium salts such as benzethonium chloride could also be used.

#### Example No.14

Different amounts of an amino-modified silicone-polyether were used to establish an acceptable range for use as a tipping lubricant.

	<u>1</u>	<u>2</u>	
	Amino-modified	1.0	10.0
	silicone surfactant (g)		
35	Water (g)	99.0	90.0
	Tip Quality	Good	Good

5

		Control		Test	
		<u>Tip Adhesion</u>	<u>Catheter Drag</u>	<u>Tip Adhesion</u>	<u>Catheter Drag</u>
		(lbs)	(g)	(lbs)	(g)
10	Aging				
	0 Week	0.19	3.0	0.12	4.8
	1 Week, 90°C	0.10	4.0	0.11	3.6
	2 Week, 90°C	0.16	3.9	0.12	4.7

From the above data, it is clear that the lubricant of this invention stabilizes tip adhesion  
 15 and the catheter lubricant lubricates the catheter adequately. Furthermore, the  
 properties of the lubricant of this invention are comparable to the control products.

#### Example No. 12

20 Other silicone surfactant combinations have been tried. For example, a  
 polyalkylene oxide-modified polydimethylsiloxane block copolymer known as Silwet  
 L7230 which is similar to Silwet L7001 was used in combination with an amino-modified  
 silicone-polyether copolymer known as Silwet Y12593.

25

		<u>Waterborne Lubricant</u>	<u>Silicone 1 MM cstk</u>
	<u>Ingredients</u>		
	Silwet L7230 (%)	4.50	-
	Silwet Y12593 (%)	0.50	-
30	Vitamin E (%)	0.25	-
	Water (%)	94.75	-
	Silicone 1 MM cstk (%)	-	2.4
	Freon TF (%)	-	97.6

35 The catheter tubing is swaged into the catheter adapter by using stainless steel  
 wedges. When the catheter assembly is put over the needle assembly it is possible  
 that the stainless steel wedge and the stainless steel needle may rub. Thus there could  
 exist a high resistive force between the two metal surfaces. This becomes even more  
 prominent when the catheter is pushed off of the needle depending upon the angle  
 40 between the catheter tubing and the needle. A test was devised in which the catheter

- 5 similarly produced using the control silicone fluid solutions. These products had good quality tips. They were packaged and sterilized by ethylene oxide and aged at 90°C for two weeks to compare their performance.

**Tipping Lubricant Composition**

<b>Ingredients</b>	<b>Composition No. 1 Percent W/W</b>	<b>Composition No. 2 Percent W/W</b>
Masil 1066C	3.50	-
Masil 29	-	3.50
Lecithin	1.00	1.00
Cosmocil	50 ppm	50 ppm
Water	95.5	95.5

**Needle Lubricant Compositions**

<b>Ingredients</b>	<b>Composition No. 1 Percent W/W</b>	<b>Composition No. 2 Percent W/W</b>
Masil 1066C	3.50	-
Masil 29	-	3.50
Vitamin E	0.25	0.25
Cosmocil	50 ppm	50 ppm
Water	96.25	96.25

**Catheter Lubricant Compositions**

<b>Ingredients</b>	<b>Composition No. 1 Percent W/W</b>	<b>Composition No. 2 Percent W/W</b>
Masil 1066C	5.5	-
Masil 29	-	5.5
Vitamin E	0.25	0.25
Cosmocil	50 ppm	50 ppm
Water	94.25	94.25

## Example No.15

The two part lubrication solution of this invention was compared to a three part lubrication solution. The following experiment compared a 4% silicone surfactant solution, and a 4% amino-modified silicone surfactant with a three part lubrication solution. This three part lubrication solution used aminopropyl terminated polydimethyl siloxane with a viscosity of 2000 centistokes as the lubricant. Such a product is available from Petrarch Systems under the PS513 trade name. 3% of the PS513 product was mixed with 1% surfactant and 96% water. The surfactant used was Tergitol S and is available from Union Carbide. The tipping results are provided below:

Parameters	PS513 Emulsion	Silicone Surfactant	Amino-modified Silicone Surfactant
Tip sticking (out of 60)	0	0	0
Tipped catheter length (in.)	2.0019(0.0010)*	2.0024(0.0014)**	1.9889(0.0008)*
Target length	2.0000	2.0000	2.0000

Note: ( ) = standard deviation  
 \* = sample size 10  
 \*\* = sample size 8

The lubrication results are provided below:

	Tip Adhesion (lbs)	Catheter Drag (g)
PS513 Emulsion	0.21 (0.02)	3.8 (0.4)
Silicone Surfactant	0.48 (0.05)	5.8 (1.0)
Amino-modified Silicone Surfactant	0.23 (0.06)	3.5 (0.3)

Note: sample size = 10  
 ( ) = standard deviation

## Example No. 16

Two lubricant solutions for use in the tipping process, lubricating the needles, and lubricating the catheters were tested simultaneously. Two thousand catheter and introducer needle assemblies were produced using one lubricant solution containing Masil 1066C, and another two thousand catheters were produced using another lubricant solution containing Masil 29. The machines ran smoothly. No significant problems surfaced, and the products produced were comparable to the control products

- 5 Performance test results before and after aging the products at 90°C for two weeks indicate that the performance of Group No. 1 and Group No. 2 products are comparable to the control.

Thus, it is seen that a new lubrication system is provided that is inexpensive, easy to control, safe, non-toxic and "environmentally friendly."

10

15

20

25

30

**Lubricant Composition for Control Products, Percent**

Ingredients	Catheter Tipping Lubricant	Needle Lubricant	Catheter Lubricant
DC360 Silicone 12.50 CSTK	-	-	2.5
IMM CSTK Silicone	-	3.0	-
Amino Modified Polydimethylsiloxane PS-513	0.65	-	-
Hydrochlorofluorocarbon (HCFC)	99.35	97.0	97.25

Products produced using Composition No. 1 are identified as Group No. 1 Products.

- 10 Similarly, those products produced using Composition No. 2 were identified as Group No. 2 Products. The performance of these products are described in Table No. 1 and Table No. 2 below.

**Table No. 1**  
**Performance of EtO Sterilized, Non-Aged Products**

Parameters	Group No. 1	Group No. 2	Control
Catheter Tip Adhesion (lbs)	0.29 (0.03)	0.21 (0.06)	0.39 (0.11)
Latex Penetration			
Needle tip (g)	16.9 (1.9)	16.6 (1.1)	15.2 (1.7)
Catheter tip (g)	16.7 (2.0)	17.6 (2.6)	11.8 (0.8)
Catheter drag (g)	4.2 (0.5)	5.2 (0.9)	3.7 (0.7)
Needle/Adapter Pull (lbs)	46.9 (3.8)	50.1 (5.0)	50.0 (4.3)
Catheter/Adapter Pull (lbs)	2.28 (0.21)	2.45 (0.24)	2.35 (0.19)

**Table No. 2**  
**Performance of EtO Sterilized, 2 Weeks Aged at 90°C Products**

Parameters	Group No. 1	Group No. 2	Control
Catheter Tip Adhesion (lbs)	0.89 (0.29)	0.26 (0.09)	1.05 (0.11)
Latex Penetration			
Needle tip (g)	21.4 (2.6)	18.1 (1.9)	13.6 (1.4)
Catheter tip (g)	17.4 (4.8)	18.8 (3.4)	12.1 (1.0)
Catheter drag (g)	3.2 (0.6)	3.7 (0.9)	3.0 (0.3)
Needle/Adapter Pull (lbs)	45.1 (3.0)	47.9 (4.3)	41.9 (4.4)
Catheter/Adapter Pull (lbs)	2.43 (0.14)	2.41 (0.12)	2.47 (0.18)

- 5           10.    The medical device of Claim 9 wherein the silicone surfactant also  
includes an amino-modified silicone-polyether copolymer.
11.    The medical device of Claim 8 wherein the silicone surfactant is an amino-  
10   modified silicone-polyether copolymer.
12.    The medical device of Claim 7 wherein the silicone based water soluble  
lubricant is an alkylene oxide-modified silicone glycol.
13.    A method for forming a shaped tip on a catheter comprising:  
15   applying a silicone based water soluble lubricant to a tip of an untipped catheter  
tubing;  
          placing the untipped catheter tubing on a mandrel;  
          heating the untipped catheter tubing to soften the tip of the tubing;  
          engaging the mandrel with a die to form the catheter tip; and  
20   removing the shaped catheter tubing from the die and mandrel.
14.    The method of Claim 13 wherein the silicone based water soluble  
lubricant is applied to a tip of an untipped catheter tubing by applying a lubricant  
25   solution of the silicone based water soluble lubricant and water to the tip.
15.    The method of Claim 14 wherein the silicone based water soluble  
lubricant is a silicone surfactant.
16.    The method of Claim 14 wherein the silicone based water soluble  
30   lubricant is an alkylene oxide modified silicone glycol.

5 We claim:

1. A lubrication solution for a medical device, comprising:  
a silicone based water soluble lubricant; and  
water.

10 2. The lubrication solution of Claim 1, wherein the silicone based water  
soluble lubricant is a silicone surfactant lubricant.

15 3. The lubrication solution of Claim 2 wherein the silicone surfactant lubricant  
is a block copolymer polyalkylene oxide-modified polydimethylsiloxane.

4. The lubrication solution of Claim 2 wherein the silicone surfactant lubricant  
is an amino-modified silicone polyether copolymer.

20 5. The lubrication solution of Claim 3 wherein the silicone surfactant lubricant  
also includes an amino-modified silicone polyether copolymer.

6. The lubrication solution of Claim 1 wherein the silicone based water  
soluble lubricant is an alkylene oxide modified silicone glycol.

25 7. A medical device having a surface with a silicone based water soluble  
lubricant deposited thereon.

30 8. The medical device of Claim 7 wherein the silicone based water soluble  
lubricant is a silicone surfactant.

9. The medical device of Claim 8 wherein the silicone surfactant is a block  
copolymer polyalkylene oxide-modified polydimethylsiloxane.

# INTERNATIONAL SEARCH REPORT

International Application No.

PCT/US 99/05846

**A. CLASSIFICATION OF SUBJECT MATTER**  
IPC 6 A61L29/00

According to International Patent Classification (IPC) or to both national classification and IPC

**B. FIELDS SEARCHED**

Minimum documentation searched (classification system followed by classification symbols)

IPC 6 A61L

Documentation searched other than minimum documentation to the extent that such documents are included in the fields searched

Electronic data base consulted during the international search (name of data base and, where practical, search terms used)

**C. DOCUMENTS CONSIDERED TO BE RELEVANT**

Category	Citation of document, with indication, where appropriate, of the relevant passages	Relevant to claim No.
X	EP 0 698 463 A (BECTON DICKINSON CO) 28 February 1996 (1996-02-28) claims	1-16
X	EP 0 709 103 A (BECTON DICKINSON CO) 1 May 1996 (1996-05-01) cited in the application abstract claims	1-12

☐ Further documents are listed in the continuation of box C.

☒ Patent family members are listed in annex.

**\* Special categories of cited documents :**

"A" document defining the general state of the art which is not considered to be of particular relevance

"E" earlier document but published on or after the international filing date

"L" document which may throw doubts on priority claim(s) or which is cited to establish the publication date of another citation or other special reason (as specified)

"O" document referring to an oral disclosure, use, exhibition or other means

"P" document published prior to the international filing date but later than the priority date claimed

"T" later document published after the international filing date or priority date and not in conflict with the application but cited to understand the principle or theory underlying the invention

"X" document of particular relevance; the claimed invention cannot be considered novel or cannot be considered to involve an inventive step when the document is taken alone

"Y" document of particular relevance; the claimed invention cannot be considered to involve an inventive step when the document is combined with one or more other such documents, such combination being obvious to a person skilled in the art.

"&" document member of the same patent family

Date of the actual completion of the international search

17 August 1999

Date of mailing of the international search report

26/08/1999


Name and mailing address of the ISA

European Patent Office, P.B. 5818 Patentlaan 2  
NL - 2280 HV Rijswijk  
Tel. (+31-70) 340-2040, Tx. 31 651 epo nt,  
Fax: (+31-70) 340-3016

Authorized officer

Muñoz, M

FIG-1

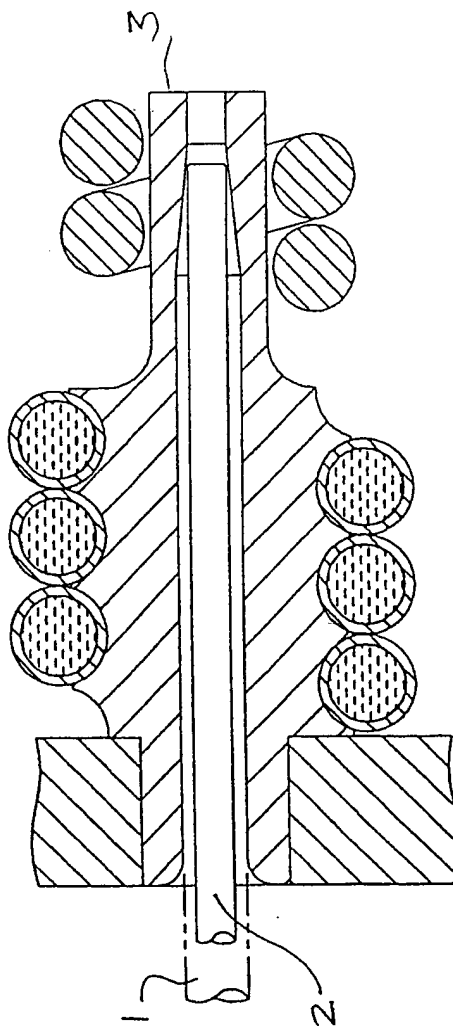

# INTERNATIONAL SEARCH REPORT

Information on patent family members

International Application No

PCT/US 99/05846

Patent document cited in search report	Publication date	Patent family member(s)	Publication date
EP 0698463 A	28-02-1996	CA 2154137 A	23-02-1996
		JP 2838989 B	16-12-1998
		JP 8066470 A	12-03-1996
		US 5589120 A	31-12-1996
		US 5688747 A	18-11-1997
EP 0709103 A	01-05-1996	CA 2154112 A	23-02-1996
		JP 8066471 A	12-03-1996
		US 5653695 A	05-08-1997

THIS PAGE BLANK (USPTO)

EPM TC 3700  
FINAL SEARCH DATE 2-11-03  
DELIVER TO GOVT DATE 2-12-03